

Endostar EASYdam Kit



Endostar EASYdam Kit

A specially selected set of tools used for placing the rubber dam.

- > Only the necessary and carefully selected elements (most popular clamps).
- > The only set on the market with the most commonly used metal frame.
- > German quality and durability.
- > Easy storage - convenient and aesthetic plastic packaging.

The kit contains:

- > **Endostar EASYdam Clamp**, rubber dam matte clamps: #9, #2, #2A, 12A, 13A 1 pc each
- > **Endostar EASYdam Forceps**
rubber dam forceps, 1 pc
- > **Endostar EASYdam Punch**
rubber dam punch, 1 pc
- > **Endostar EASYdam**
latex sheets, medium, blue, for adults, odorless, 36 pcs
- > **Endostar EASYdam Template**
rubber dam plastic template, 1 pc
- > **Endostar EASYdam Frame**
metal dam frame for adults, 1 pc
- > **Endostar EASYfix**
Ø 1,5 mm, yellow, small, length: 4 m, 1 pc




Poldent Co. Ltd.
Dzika 2 Street, 00-194 Warsaw, Poland
Phone: +48 22 351 76 50 - 51
E-mail: poldent@poldent.pl
www.poldent.pl

Product information

EASYdam - latex rubber dam for adults			
EEDLM36G	152 x 152 mm, medium, odorless, green, 36 pcs		
EEDLH36G	152 x 152 mm, thick, odorless, green, 36 pcs		
EEDLT36G	152 x 152 mm, thin, odorless, green, 36 pcs		
EEDLM36B	152 x 152 mm, medium, odorless, blue, 36 pcs		
EASYdam NON LATEX - nonlatex rubber dam for adults			
EEDNL30B	150 x 150 mm, medium, odorless, blue, 30 pcs		
Endostar EASYdam Clamp - matte clamps with wings			
EEDC9	#9	EEDC14A	#14A
EEDC2	#2	EEDC7	#7
EEDC2A	#2A	EEDC3	#3
EEDC12A	#12A	EEDC8A	#8A
EEDC13A	#13A		
Endostar EASYdam Frame			
EEDFR	Metal frame for adults, 1 pc		
Endostar EASYdam Punch			
EEDP	Punch, 1 pc		
Endostar EASYdam Forceps			
EEDFO	Forceps, 1 pc		
Endostar EASYdam Template			
EEDT	Plastic template, 1 pc		
Endostar EASYdam Napkins			
EEDN50	Napkins, 50 pcs		
Endostar EASYfix - rubber cord for stabilization			
EEDEFS	Ø 1,5 mm, yellow, small, length 4 m, 1 pc		
EEDEFL	Ø 2,0 mm, orange, large, length 4 m, 1 pc		

Endostar EASYdam

Make it EASY

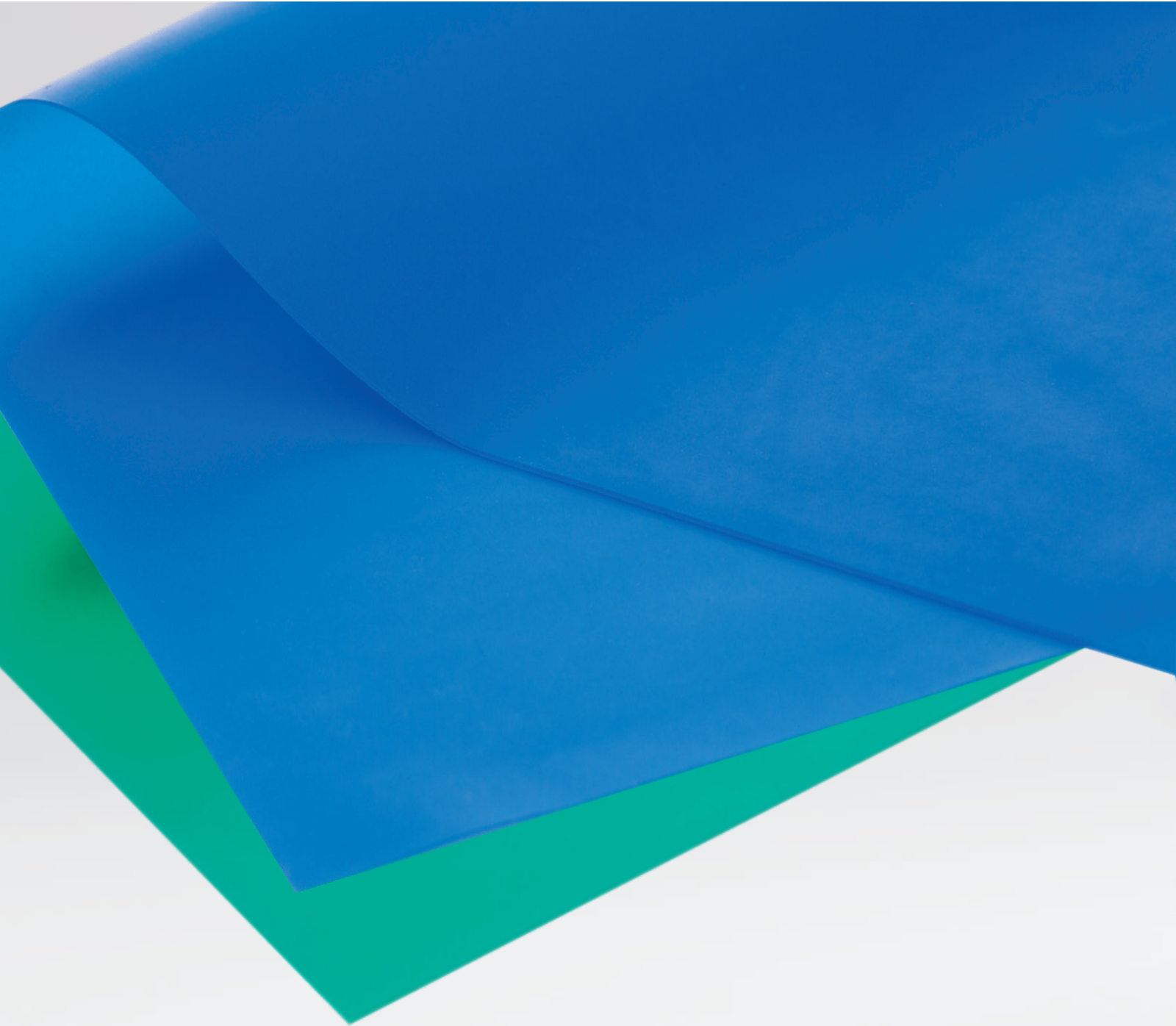

www.poldent.pl



Endostar EASYdam

Endostar EASYdam

The extremely durable and flexible powder-free rubber dam provides comfort for the patient. Very flexible and with high resistance to tearing. Available in two colors: blue and green and three thickness options: thin, medium, thick. Odorless. Very easy to apply.


Endostar EASYdam Non-latex

Latex-free, powderless rubber dam with increased flexibility compared to a traditional latex rubber dam. Extremely easy to apply. Much more flexible than traditional latex rubber dam. Odorless and tasteless. Closely adheres to the tooth. Non-allergenic.



Endostar EASYdam Clamps

High quality matte rubber dam clamps with wings. The matte finish guarantees a comfortable work environment without reflections. This is particularly useful when working with a microscope. The high-quality steel provides high wear resistance and fatigue strength.

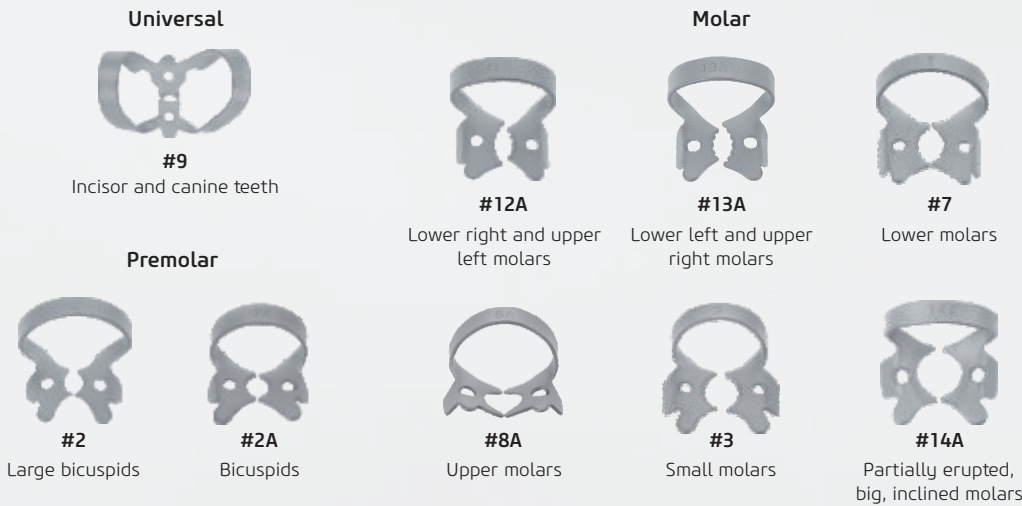

EASYdam accessories



Endostar EASYdam Frame

The adult-sized metal frame firmly holds the rubber dam in place. U-shaped frame. Autoclavable at 134°C.



Endostar EASYdam Punch

The Ainsworth punch provides precise punching power for each rubber dam thickness option in 5 different sizes. Suitable for frontal and posterior teeth. The highest quality of stainless steel instruments guarantees precise puncturing.



Endostar EASYdam Forceps

Comfortable rubber dam forceps allow quick and safe transportation, application and removal of clamps. Made from stainless steel. Matte finish.




Endostar EASYdam fix

The Easydam Fix provides support or can be an alternative for traditional clamps. Single-use. In a practical container. Easydam Fix occurs in length 4 m and 2 thickness options: 2mm (orange) and 1.5 mm (yellow).



Endostar EASYdam Template


The Endostar EASYdam template helps accurately determine the position of teeth on the rubber dam sheet. Very helpful in proper punching. Plastic and durable and easy to use.

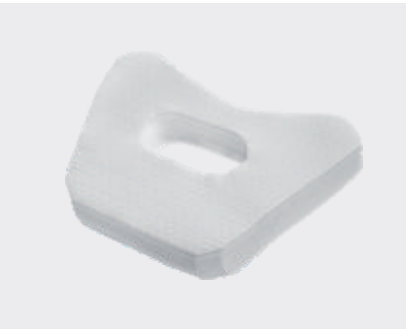


Endostar EASYdam Napkins

Rubber dam napkins for maximum patient comfort. Napkins absorb moisture and protect the patient's skin.
> Protects delicate skin.
> Absorbs water, saliva and sweat.
> Provides a dry contact surface.
> 50 pieces.

Step by step

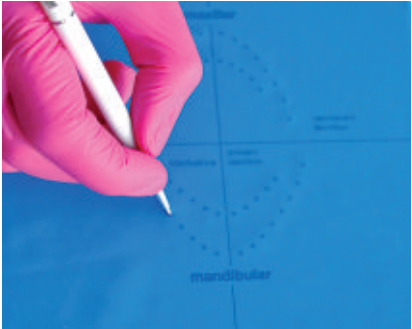

The Endostar EASYdam Templates provides clear and simple markings.

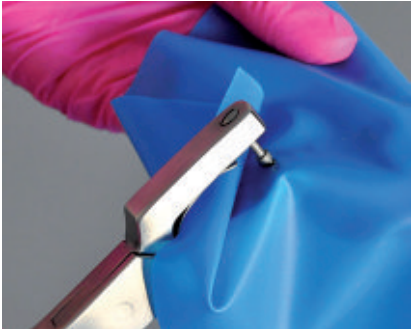

The Endostar EASYdam Punch for perforating the rubber dam. Precise punching without tearing.




Place the clamp with the rubber dam on the tooth and use the frame to secure the borders of the rubber dam.



Insert the Endostar EASYdam Clamp into the rubber dam with the Endostar EASYdam Forceps.



Always ensure the clamp with dental floss before placing it in the mouth of the patient.

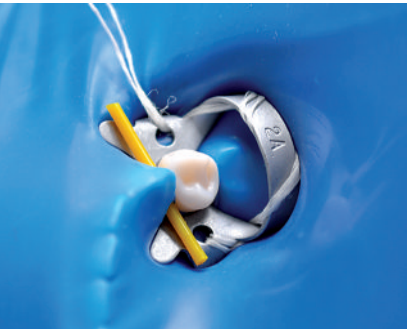

Endostar EASYfix is a supplement or alternative to traditional clamps.