

Complementary products

To improve endo treatments



Endostar ENDObox without/with instruments

ENDObox - special box for holding endodontic instruments and for their sterilization.

Endostar ENDObox without/with instruments

ENDObox is used for holding endodontic instruments and for their sterilization. A two piece box (bottom and top part) is made out of anodized aluminium and has a movable endodontic insert.

ENDObox with instruments consists of 28 endodontic instruments for root canal treatment:

- > 9 rotary files - Endostar E3 Rotary System (Basic, Big Apical, Small Apical)

- > 12 K-files with endostops, length 25 mm and sizes: 06 rose, 08 grey, 10 violet, 15 white, 20 yellow, 25 red
- > 6 NiTi K-files with endostops, length 25 and sizes: 25 red, 30 blue, 35 green, 40 black, 45 white
- > 1 Endostar Canal locator, length 18 and size: 10 violet

Sterilization in steam autoclave up to 134°C.

Please ensure any solution used is checked for compatibility with aluminium.



File type:

Endostar ENDObox

Endostar K-files	06 x 2 / 08 / 10 x 3 / 15 x 3 / 20 x 2 / 25
Endostar NiTi K-files	25 / 30 x 2 / 35 / 40 / 45
Endostar E3 Small	20/06, 25/04, 20/04
Endostar E3 Big	35/04, 40/04, 45/04
Endostar E3 Basic	30/08, 25/06, 30/04
Endostar Canal Locator	10

Endostar Paper/Gutta Percha Points

Complementary products and comprehensive solutions. Endostar points with the length marking. Calibration saves time of every endo treatment.

Endostar Paper Points

Paper points for drying the root canal with a very high absorption capacity.

- > Used for drying root canal.
- > ISO coded.
- > Highly absorbent.
- > Firm and elastic.
- > Available in 12 sizes (15, 20, ..., 80) and/or assorted 15-40, 45-80, packed in compartment boxes, with standard 2% taper, 200 pcs/box or with taper 4% and 6%, 100 pcs/box.

Endostar Gutta Percha Points

For root canal obturation.
ISO coded.

- > Ideal for vertical and lateral root canal obturation.
- > Biocompatible and hygienic.
- > With high dimensional stability.
- > Firm and elastic.
- > Available in 12 sizes (15, 20, ..., 80) and/or assorted 15-40, 45-80, packed in compartment boxes, with standard 2% taper, 120 pcs/box or with taper 4% and 6%, 60 pcs/box.
- > Calibrated, with the length marking on the point.



Points type:	Paper Points	Gutta Percha Points
Size / ISO:	15-80	15-80
Taper:	2% / 4% / 6%	2% / 4% / 6%
Packing:	Single size or Assorted	Single size or Assorted
Quantity in package:	200 pcs / 100 pcs	120 pcs / 60 pcs

Essential supplement

Products for better and ergonomic treatment.



Endostar ENDOstand

Hinged stand with lid is used for holding endodontic files and reamers. Stand with a ruler to measure stop positions easily from either direction.

Sterilization in steam autoclave up to 134°C.




Endostar ENDOrepair MTA

Endostar ENDOrepair MTA is a powder consisting of fine, hydrophilic particles that set in the presence of water. Hydration of the powder creates a colloidal gel that solidifies to form a strong impermeable barrier that fully cures over a four - week period.

Benefits

- > Excellent biocompatibility.
- > Hydrophilic.
- > Water - Based Chemistry.
- > Calcium ions release.
- > Tooth coloured formula.
- > High radiopacity.

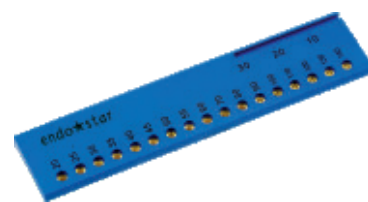

Endostar ENDOcalibrator

Gauge is used to customize accessory sized gutta percha to ISO sizes, and to measure files, gutta percha points and absorbent points.

Sterilization in steam autoclave up to 132°C.

Poldent Co. Ltd.
Dzika 2 Street
00-194 Warsaw, Poland
Mob.: +48 603 630 090
Ph.: +48 22 351 76 61
e-mail: endostar@endostar.eu
www.endostar.eu

endo★star