

Complementary products

To improve endo treatments



ENDO-irrigation

Effective, easy and safe irrigation of root canals.

Endostar ENDONEEDLES


- > Single-use disposable endo-irrigation needles with single or double vent.
- > Sterile.
- > Pyrogen-free.
- > Application: rinsing root canals and also for rinsing gingival pockets.
- > The needle is made of the highest quality material, it is characterized by thin walls and high flexibility allowing for adaptation of the needle shape to the root canal.
- > Available sizes:
 - endo-irrigation needles with single vent, size 0,3 x 25 (30G x 1")
 - endo-irrigation needles with single vent, size 0,5 x 25 (25G x 1")

- endo-irrigation needles with double vent, size 0,3 x 25 (30G x 1")

Endostar ENDOSYRINGE

Small, colored 5 ml irrigation syringes.

- > Organization and identification
 - color marking has a decisive impact on avoiding a mistake in changing syringes with rinsing solutions.
- > Increasing the safety of rinsing (available in 4 colors: transparent, blue, red, yellow).
- > Color coding for convenient work with the microscope.
- > Available in 5 ml capacity.
- > Standard luer-lock.

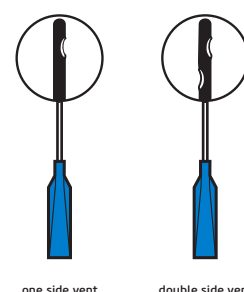

Points type:	Endostar ENDONEEDLES	Endostar ENDOSYRINGE
Size:	0,3 x 25 (30G x 1") / 0,5 x 25 (25G x 1") / 0,3 x 25 (30G x 1")	5 ml
Packing:	Box	Box
Quantity in package:	100 pcs	60 pcs

Endostar Paper/Gutta Percha Points

Complementary products and comprehensive solutions. Endostar points with the length marking. Calibration saves time of every endo treatment.

Endostar Paper Points

Paper points for drying the root canal with a very high absorption capacity.

- > Used for drying root canal.
- > ISO coded.
- > Highly absorbent.
- > Firm and elastic.
- > Available in 12 sizes (15, 20, ..., 80) and/or assorted 15-40, 45-80, packed in compartment boxes, with standard 2% taper, 200 pcs/box or with taper 4% and 6%, 100 pcs/box.

Endostar Gutta Percha Points

For root canal obturation.
ISO coded.


- > Ideal for vertical and lateral root canal obturation.
- > Biocompatible and hygienic.
- > With high dimensional stability.
- > Firm and elastic.
- > Available in 12 sizes (15, 20, ..., 80) and/or assorted 15-40, 45-80, packed in compartment boxes, with standard 2% taper, 120 pcs/box or with taper 4% and 6%, 60 pcs/box.
- > Calibrated, with the length marking on the point.



CE
0197

Points type:	Paper Points	Gutta Percha Points
Size / ISO:	15-80	15-80
Taper:	2% / 4% / 6%	2% / 4% / 6%
Packing:	Single size or Assorted	Single size or Assorted
Quantity in package:	200 pcs / 100 pcs	120 pcs / 60 pcs

Essential supplements

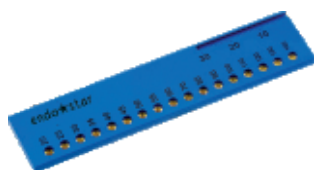


Endostar ENDocalibrator

Gauge is used to customize accessory sized gutta percha to ISO sizes, and to measure files, gutta percha points and absorbent points. Sterilization in steam autoclave 134°C.



Endostar ENDObox

Endodontic container designed for storing and sterilizing instruments used in root canal treatment. It has 30 slots for rotary or hand files. Sterilization in steam autoclave 134°C.

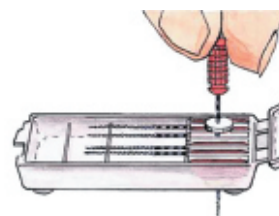

Endostar Endostops

Round, white silicon rubber, stoppers, used for marking the working length on a file or reamer. Easily visible on RTG. Sterilization in steam autoclave 134°C.



Endostar ENDostand

Hinged stand with lid holds endodontic files and reamers. Stand with a ruler to measure stop positions easily from either direction. Sterilization in steam autoclave 134°C.



Endostar ENDosizer

A special instrument designed for precise measurement of the working length of files and the length of gutta-percha and paper points. It is made of extremely durable polycarbonate and consists of a recess for endostop and a separate canal for measuring and cutting off points. Sterilization in steam autoclave 134°C.

Ver. 2, December 2023

Poldent Co. Ltd.
Dzika 2 Street
00-194 Warsaw, Poland
Ph.: +48 22 351 76 61
e-mail: endostar@endostar.eu
www.endostar.eu

endo★star